

Stand Out with Your Excellent Cosmetic Advice

Multi Skin Test Center[®] MC 1000 & Software

moisture, oil, melanin, pH, elasticity & many more




Allows operating all C+K probes at the same time according to individual demands due to modular C+K probe technology.

C+K-measuring probes for following parameters can be connected: moisture, oil, pH, pigmentation & redness, transepidermal waterloss (TEWL – the barrier function of the skin), skin temperature and elasticity.

Additionally the cameras Visioscope[®] PC35 or/and Visiopor[®] PP 34N can be connected.



The **Multi Skin Test Center[®] MC750** offers the features like the MC 1000 except for the camera functions.

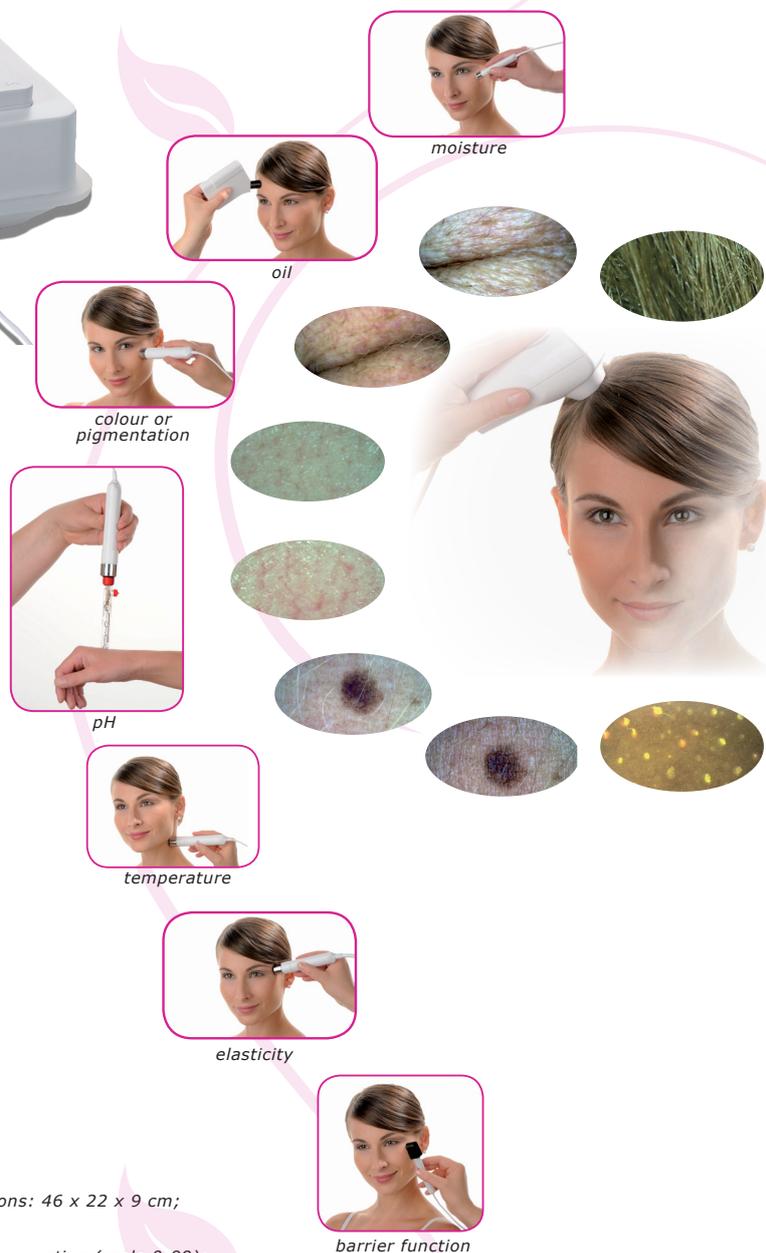

Please ask for our special transportation accessories.

Technical Data:

External power supply: Input: 100-240 V, 50-60 Hz, Output: 12 V, 4 A, Dimensions: 46 x 22 x 9 cm;
Weight without probes: approx. 1.8 kg, Interface: USB 2.0

Measuring principles:

Moisture: capacitance (scale 0-99); Sebum: photometrical (scale 0-99); Elasticity: suction (scale 0-99);
Melanin/Erythema: absorption/reflection (scale 0-99); pH: glass electrode (pH with one decimal);
TEWL-Index: open chamber water flux measurement (scale 0-20)

Additional Visioscope[®] PC35 or Visiopor[®] PP 34N

Devices must be operated by an IBM-compatible computer with Windows[®] XP or higher and a USB-port.
Technical changes may be made without prior notice.







www.courage-khazaka.de

Cologne Germany

Stand Out with Your Excellent Cosmetic Advice

Multi Skin Test Center® MC 1000 & Software

moisture, oil, melanin, pH, elasticity and many more

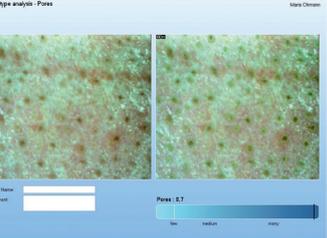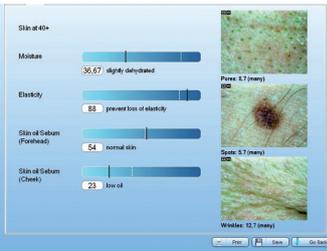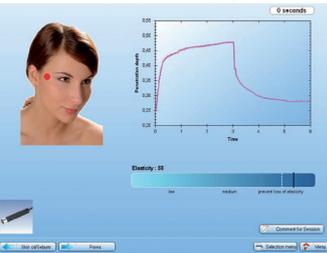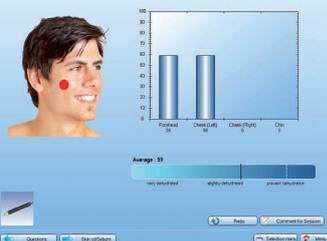

- ❖ **CSI-Software**
Very user-friendly and appealing software "Complete Skin Investigation" simplifying work with the system
- ❖ Individual product recommendations for skin care, hair care and sun screen products which can be printed out for the customer
- ❖ Easy increase of sales by offering a special service to distinguish from competitors and to secure customers retention

The software offers 6 modules:

- ❖ **Visioface® 1000D** (see page 15/16) programme is integrated
- ❖ Questionnaire combined with moisture, sebum and elasticity measurements as well as the camera measurements (optional) calculates **skin type** and recommends a set of skin care products (from cleansing to day and night care, additional care and food supplements) for women or men
- ❖ **Imaging** with Visioscope® PC35 and calculation of pores, spots and wrinkles. Easy comparison of customer images side by side. The Visiopor® PP34N shows the porphyrines from the acne bacteria.

- ❖ Sebum production and desquamation can be analyzed with the special foils Sebu-fix® and Corneofix®.
- ❖ **Hair analysis** (measurement of sebum, colour, imaging for dandruff, hair thickness) with questionnaire and product recommendation
- ❖ Answering questions and measurement of melanin content allows determination of most suitable **sun screen** and aftersun products, taking into account where, when and how long somebody is going to be exposed to the sun including the respective UV-index
- ❖ **Single measurements** (moisture, sebum, elasticity, pigmentation/erythema, pH, temperature and TEWL) in any order and on any body site to be saved and printed out

- ❖ User can easily enter the products to be recommended in each programme according to the age of the person, the measurement results and answers to the questions
- ❖ The appearance of the software can be easily changed by altering the name of the programme, exchange colours or the start screen picture, add a logo or even create a new questionnaire.

